

# ***Gynaecological Essentials***

*Incorporating  
Uterine Aspiration Essentials*

**purple surgical**





# Secondary Trocars

The Purple Surgical Secondary Trocar range consists of single-use 5mm ports in a classic and proven design. These low profile trocars allow controlled entry to the abdomen under visual supervision with no tenting of the peritoneum. Atraumatic retention threads on the cannula ensure minimum protrusion into the cavity and maximum working space in the pelvis. All versions carry a closed luer connector for insufflation filter and tubing set connection as required.

Features and benefits include:

- Smooth obturator to cannula transition and honed stainless steel tips for precise, controlled entry
- Retention threads ensure minimal cavity protrusion and smooth instrument exchange
- One piece silicone seal guarantees minimal leaks with or without instruments
- Available in standard or extended length, with conical or pyramidal tips



*Pyramidal Tip*

*Conical Tip*

Product	Code	Description	Pack Qty
Secondary Trocar	PS9800	5mm x 55mm Pyramidal Tip Trocar with Gas Port	5 / Box
Secondary Trocar	PS9801	5mm x 55mm Conical Tip Trocar with Gas Port	5 / Box
Secondary Trocar	PS9840	5mm x 70mm Pyramidal Tip Trocar with Gas Port	5 / Box
Secondary Trocar	PS9841	5mm x 70mm Conical Tip Trocar with Gas Port	5 / Box
Secondary Trocar Kit	PS9820	5mm x 55mm Pyramidal Tip Trocar with Gas Port + Additional Cannula	5 / Box
Secondary Trocar Kit	PS9821	5mm x 55mm Conical Tip Trocar with Gas Port + Additional Cannula	5 / Box
Secondary Trocar Kit	PS9842	5mm x 70mm Pyramidal Tip Trocar with Gas Port + Additional Cannula	5 / Box
Secondary Trocar Kit	PS9847	5mm x 70mm Conical Tip Trocar with Gas Port + Additional Cannula	5 / Box



## Spackman Cannula

The Purple Surgical Spackman Cannula allows for the intra-operative manipulation of both the anteverted and retroverted uterus as well as facilitating quick, simple, effective hysterosalpingogram for tubal patency studies.

Features and benefits include:

- Cost effective single-use device eliminates the need for cleaning and sterilisation
- Clearly marked with universal single-use logo to ensure disposal after use
- High quality, medical grade stainless steel
- Weighted handle provides effective 'hands-free' anteversion
- Anatomical cervical cone ensures effective, leak free seal
- Convenient, adjustable tenaculum attachment
- The atraumatic tip provides reduced uterine perforation risk




Product	Code	Description	Pack Qty
Spackman Cannula	PS3112	Spackman Uterine Manipulator & Injector	25 / Box

## Purple Uterine Manipulator & Injector (PUMI)

The Purple Surgical PUMI allows for the intra-operative manipulation of both the anteverted and retroverted uterus as well as facilitating quick, simple, effective hysterosalpingogram for tubal patency studies.

Features and benefits include:

- Anatomic balloon cuff ensures effective, leak free seal whilst eliminating uterine spasm and cervical trauma due to tenaculum use
- Rigid, curved handle for both the anteverted and retroverted uterus
- Convenient alternative to the difficult to clean and sterilize, reusable Spackman style cannula
- Adjustable handle with centimetre markings enable precise uterine depth setting that can be adapted to fit the patient's anatomy
- Atraumatic tip provides reduced uterine perforation risk
- Minimal dilation is required due to the 4.5mm cannula diameter

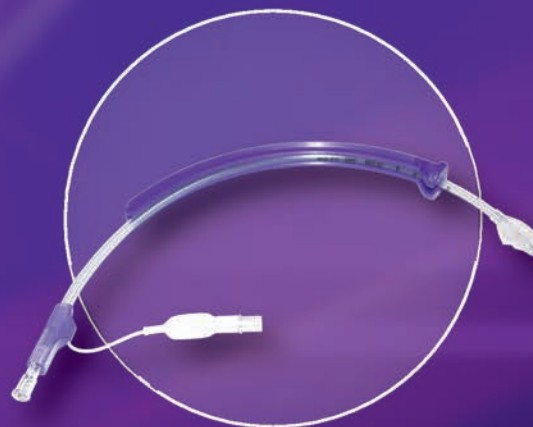

Product	Code	Description	Pack Qty
PUMI	PS3111	Purple Uterine Manipulator & Injector	10 / Box




## Rectal and Vaginal Visualisation Probes

The Purple Surgical Rectal and Vaginal Visualisation Probes provide effective surgeon guidance. Complex pelvic adhesiolysis, endometriosis and laparoscopic hysterectomy procedures require the precise localisation and visualisation of both vaginal and rectal structures.

Features and benefits include:

- Rounded, atraumatic tips for precise location of rectal and vaginal structures
- Ergonomic design makes for easy introduction and use
- Integrated visual and tactile position indicators
- Robust, simple to maintain, reusable construction

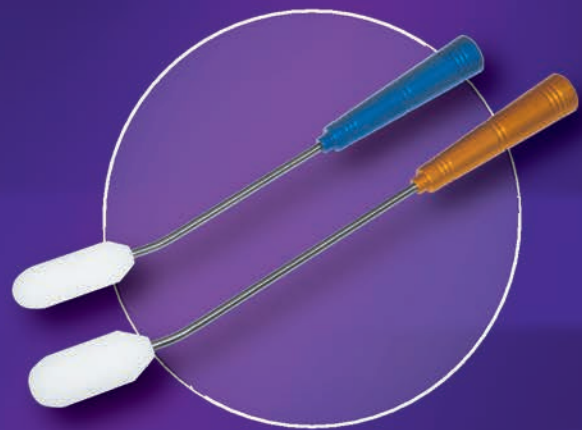


Product	Code	Description	Colour	Pack Qty
Visualisation Probe	PS9590	Reusable Vaginal Probe	Blue	1 / Box
Visualisation Probe	PS9595	Reusable Rectal Probe	Orange	1 / Box

## Linear wire Cone Biopsy Electrodes (LCBE)

The Purple Surgical single-use Linear wire Cone Biopsy Electrode range comprises of 5 different diameters and depths for easy, controlled biopsies providing true cone samples with improved margin quality for better analysis.

Features and benefits include:

- Extended tip locates in the cervical os for optimal stability
- Single-use for convenience, hygiene and reliability
- Linear wire provides stable current characteristics and better quality samples
- Curved and inclined guide arm ensures depth control and snag-free clockwise rotation
- Adhesive sample orientation label included
- Strong tungsten cutting wire minimises incidence of wire breakage

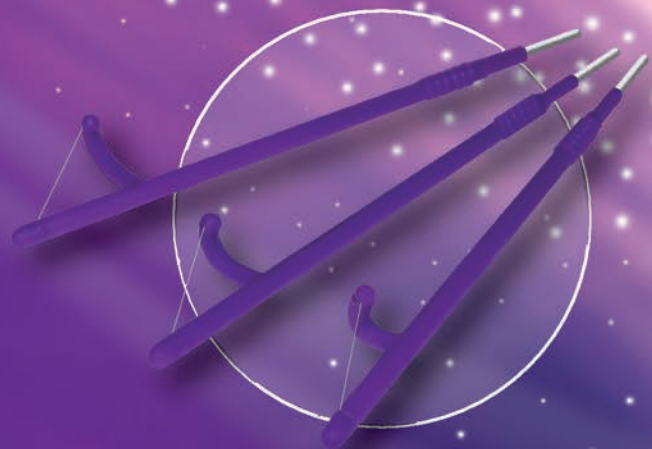

Product	Code	Description	Pack Qty
Linear wire Cone Biopsy Electrodes (LCBE)	PS9150	Small (9mm Radius x 12mm Depth)	5 / Box
Linear wire Cone Biopsy Electrodes (LCBE)	PS9151	Medium (12mm Radius x 15mm Depth)	5 / Box
Linear wire Cone Biopsy Electrodes (LCBE)	PS9152	Large (15mm Radius x 18mm Depth)	5 / Box
Linear wire Cone Biopsy Electrodes (LCBE)	PS9154	Medium Deep (12mm Radius x 20mm Depth)	5 / Box
Linear wire Cone Biopsy Electrodes (LCBE)	PS9155	Wide Shallow (15mm Radius x 8mm Depth)	5 / Box
Linear wire Cone Biopsy Electrodes (LCBE)	PS9156	Multi Set (1x Each Size)	5 / Box




## LEEP Electrodes

Purple Surgical's LEEP Electrodes are high quality, single-use electrodes that guarantee performance, every time.

Features and benefits include:

- 9 different combinations of width and depth provide all of the most frequently used electrode dimensions
- 1 square box electrode for 'top hat' excisions
- 1 extended length ball electrode for effective haemostasis of the excision site

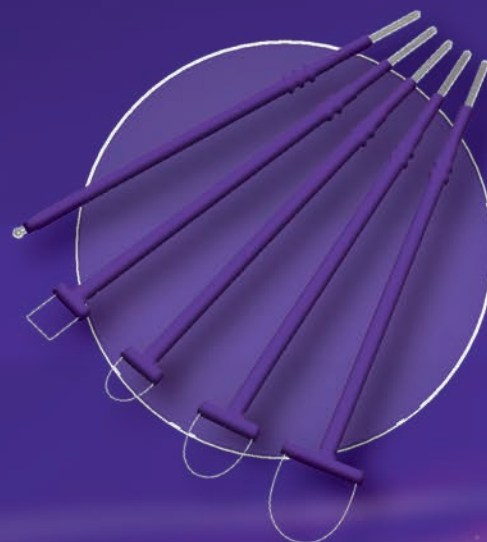


Product	Code	Description	Pack Qty
LEEP Electrode	PS9160	Small Loop Electrode with 12cm Shaft (1cm Width x 0.7cm Depth)	5 / Box
LEEP Electrode	PS9161	Small Loop Electrode with 12cm Shaft (1cm Width x 1cm Depth)	5 / Box
LEEP Electrode	PS9162	Medium Loop Electrode with 12cm Shaft (1.5cm Width x 0.5cm Depth)	5 / Box
LEEP Electrode	PS9163	Medium Loop Electrode with 12cm Shaft (1.5cm Width x 0.7cm Depth)	5 / Box
LEEP Electrode	PS9164	Medium Loop Electrode with 12cm Shaft (1.5cm Width x 1.5cm Depth)	5 / Box
LEEP Electrode	PS9165	Large Loop Electrode with 12cm Shaft (2cm Width x 0.8cm Depth)	5 / Box
LEEP Electrode	PS9166	Large Loop Electrode with 12cm Shaft (2cm Width x 1cm Depth)	5 / Box
LEEP Electrode	PS9167	Large Loop Electrode with 12cm Shaft (2cm Width x 1.5cm Depth)	5 / Box
LEEP Electrode	PS9168	Large Loop Electrode with 12cm Shaft (2cm Width x 2cm Depth)	5 / Box
LEEP Electrode	PS9169	Square Electrode with 12cm Shaft (1cm Width x 1cm Depth)	5 / Box
LEEP Electrode	PS9170	Ball Electrode with 12cm Shaft (4mm Diameter)	5 / Box

## Abrasive Tip Cleaner

The Purple Surgical Abrasive Tip Cleaner allows for the effective removal of charred debris from electrodes.

Features and benefits include:

- Large 5cm x 5cm effective area
- Extended release liner enables convenient removal with gloved hands
- Radio-opaque marker ensures patient safety through rapid detection
- Coated abrasive surface ensures no flaking of contaminated abrasive material during use

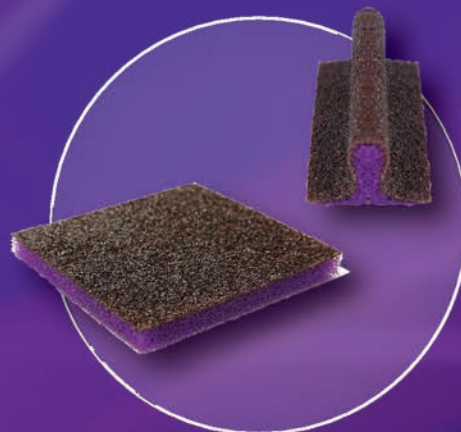

Product	Code	Description	Pack Qty
Abrasive Tip Cleaner	PS5200	5cm x 5cm	25 / Box
Abrasive Tip Cleaner	PS5204	5cm x 5cm	100 / Box

## Handswitch Operated Electrosurgical Pencil

The Purple Surgical Handswitch Operated Electrosurgical Pencil is a reliable and cost-effective single-use pencil for LEEP type electrode procedures.

Features and benefits include:

- Sterile, single-use convenience
- Comfortable, ergonomic design
- Tactile cut and coag button action
- Universal ESU fitting



Product	Code	Description	Pack Qty
Electrosurgical Pencil	PS5555	Handswitch Operated Pencil	25 / Box

## Footswitch Operated Electrosurgical Pencil

The Purple Surgical Footswitch Operated Electrosurgical Pencil is a practical and reliable single-use pencil for footswitch operation allowing easy rotation of linear wire type cone biopsy electrodes.

Features and benefits include:

- Sterile, single-use convenience
- Comfortable, ergonomic design
- Footswitch operation frees hand movement for electrode rotation
- Universal monopolar ESU connector



Product	Code	Description	Pack Qty
Electrosurgical Pencil	PS5550	Footswitch Operated Pencil	25 / Box



## Rigid Curettes – Straight

The Purple Surgical Rigid Straight Curettes are designed for safe, effective uterine aspiration.

Features and benefits include:

- Reduced perforation risk due to the smooth atraumatic tip and depth gauge markings
- Provides quick, simple aspiration port orientation by means of the linear markings



Product	Code	Description	Pack Qty
Rigid Curettes - Straight	PS6007	7mm Rigid Straight	10 / Box
Rigid Curettes - Straight	PS6008	8mm Rigid Straight	10 / Box
Rigid Curettes - Straight	PS6009	9mm Rigid Straight	10 / Box
Rigid Curettes - Straight	PS6010	10mm Rigid Straight	10 / Box
Rigid Curettes - Straight	PS6011	11mm Rigid Straight	10 / Box
Rigid Curettes - Straight	PS6012	12mm Rigid Straight	10 / Box

## Rigid Curettes – Curved

The Purple Surgical Rigid Curved Curettes are designed for safe, effective uterine aspiration.

Features and benefits include:

- Reduced perforation risk due to the smooth atraumatic tip and the depth gauge markings
- Anatomic curved design to facilitate access to the anteverted uterus
- Provides quick, simple aspiration port orientation by means of the linear markings



Product	Code	Description	Pack Qty
Rigid Curettes - Curved	PS6107	7mm Rigid Curved	10 / Box
Rigid Curettes - Curved	PS6108	8mm Rigid Curved	10 / Box
Rigid Curettes - Curved	PS6109	9mm Rigid Curved	10 / Box
Rigid Curettes - Curved	PS6110	10mm Rigid Curved	10 / Box
Rigid Curettes - Curved	PS6111	11mm Rigid Curved	10 / Box
Rigid Curettes - Curved	PS6112	12mm Rigid Curved	10 / Box



## Flexible (Karmann) Curettes

The Purple Surgical Flexible Curettes are designed for early term, outpatient aspiration with minimised perforation risk.

Features and benefits include:

- Minimised perforation risk due to the soft flexible material with an atraumatic tip and depth gauge markings
- Allows the curette to collapse when pressed against the uterine fundus, due to the two distal openings in 4-8mm sizes, further reducing perforation risk and improving patient comfort
- Offers a complete range with 4, 5, 6, 7, 8, 10 and 12mm size options
- Provides quick, simple aspiration port orientation by means of the linear markings



Product	Code	Description	Pack Qty
Flexible (Karmann) Curettes	PS6204	4mm Flexible	10 / Box
Flexible (Karmann) Curettes	PS6205	5mm Flexible	10 / Box
Flexible (Karmann) Curettes	PS6206	6mm Flexible	10 / Box
Flexible (Karmann) Curettes	PS6207	7mm Flexible	10 / Box
Flexible (Karmann) Curettes	PS6208	8mm Flexible	10 / Box
Flexible (Karmann) Curettes	PS6210	10mm Flexible	10 / Box
Flexible (Karmann) Curettes	PS6212	12mm Flexible	10 / Box

## Collection Sets

The Purple Surgical Collection Sets are thick walled, high vacuum, sterile tubing with unique swivel handle for connection to sterile curettes.

Features and benefits include:

- Thick wall tubing prevents tubing wall to collapse under high vacuum pressure
- Ensures surgeon comfort and convenience due to 360° swivel handle design
- Provides optimal control due to the incorporation of a finger slide vacuum control
- Offers comfort and convenience due to the unobtrusive low torque tubing
- Available in standard 10mm (3/8") diameter with male and female connections



Product	Code	Description	Pack Qty
Collection Sets	PS6100	2m x 10mm Tubing – Male	10 / Box
Collection Sets	PS6101	2m x 10mm Tubing – Female	10 / Box



## Specimen Cup

The Purple Surgical Specimen Cups safely and effectively separate tissue from fluid, thus facilitating a 'closed' system by eliminating the need to handle aspirated tissue.

Features and benefits include:

- Provides tissue separation by means of an integral filter
- May be used for up to 12 weeks gestation, providing a solution for the majority of uterine aspiration procedures
- Allows a quick, effective visual contents check, due to the clear lid construction
- Offers safe, leak free storage and transportation
- Allows for convenient disposal by incineration and cremation
- Allows for the simple and safe addition of formaldehyde
- Ensures hygienic transportation due to the supply of re-sealable bag
- Available in an opaque version for sensitive transportation
- Available with optional box for sensitive transportation and disposal



Product	Code	Description	Pack Qty
Specimen Cups	PS6000	Medium Clear up to 12 Weeks	10 / Box
Specimen Cups	PS6666	Medium Opaque up to 12 Weeks	10 / Box
Specimen Cups	PS6000B	Medium Clear with Box	10 / Box
Specimen Cups	PS6666B	Medium Opaque with Box	10 / Box

## SafeTerm Canister

The Purple Surgical SafeTerm Canister completes the Purple Surgical 'Closed System', by providing safe, effective fluid collection with the option of an integral tissue filter in one device.

Features and benefits include:

- Simple, safe options for fluid collection alone (when used with Purple Surgical Specimen Cups) or tissue separation with the incorporation of an integral filter
- Ideal for later stage terminations up to 16 weeks
- Offers maximum flexibility due to the incorporation of 10mm male, 10mm female, 12mm female, Specimen Cup and standard suction tubing ports
- Allows rapid, accurate fluid loss assessment due to the clear volume markings
- Allows for convenient disposal by incineration and cremation
- Conveniently attaches to existing vacuum pumps by means of the universal bracket
- Minimised fluid spillage risk with the option of pre-loaded solidifying granules



Product	Code	Description	Pack Qty
SafeTerm Canister	PS6900	900cc SafeTerm Canister	10 / Box
SafeTerm Canister	PS6901	900cc SafeTerm Canister with Filter	10 / Box
SafeTerm Canister	PS6902	900cc SafeTerm Canister and Solidifying Granules	10 / Box
SafeTerm Canister	PS6903	900cc SafeTerm Canister with Filter and Solidifying Granules	10 / Box
SafeTerm Canister	PS6180	1800cc SafeTerm Canister	10 / Box
SafeTerm Canister	PS6181	1800cc SafeTerm Canister with Filter	10 / Box
SafeTerm Canister	PS6182	1800cc SafeTerm Canister and Solidifying Granules	10 / Box
SafeTerm Canister	PS6183	1800cc SafeTerm Canister with Filter and Solidifying Granules	10 / Box





PS-GYN-LUT-005-07-15



Purple Surgical International Limited | Purple Surgical UK Limited  
2 Chestnut House, Farm Close, Shenley, Herts WD7 9AD, England  
Tel: +44 (0)1923 839333 Fax: +44 (0)1923 839444 Email: [mail@puplesurgical.com](mailto:mail@puplesurgical.com)

[www.puplesurgical.com](http://www.puplesurgical.com)